

October 12, 2014

MetAmp: a tool for Meta-Amplicon analysis User Manual

Ilya Y. Zhbannikov¹, Janet E. Williams¹, James A. Foster^{1,2,3}

³Institute for Bioinformatics and Evolutionary Studies, University of Idaho, Moscow, ID 83844-3051,

²Department of Biological Sciences, University of Idaho, Moscow, ID 83844-3051,

¹Department of Bioinformatics and Computational Biology, University of Idaho, Moscow, ID 83844-3051,

Keywords:

Bioinformatics, 16S genes, microbial studies, meta-analysis, data mining

Correspondence:

Ilya Y. Zhbannikov

Department of Bioinformatics and Computational Biology

University of Idaho

Moscow, ID 83844-30521

Email: ilyaz@uidaho.edu, foster@uidaho.edu

1 Description

MetAmp tool was developed to analyze microbial amplicon data by combining several marker regions from 16S rRNA gene. Such marker regions serve as unique identifications for species. There are nine marker regions in total in bacterial 16S rRNA gene.

MetAmp was developed by Ilya Y. Zhbannikov (ilyaz@uidaho.edu, i.zhbannikov@mail.ru), James A. Foster (foster@uidaho.edu) and Janet Williams (janetw@uidaho.edu).

2 Installation

MetAmp is easy to install. However, you must have R (www.r-project.org), Python (www.python.org) and GCC (<https://gcc.gnu.org/>).

1. Download or clone MetAmp from our GitHub repository: <http://github.com/izhbannikov/MetAmp>
2. In Terminal: `cd` to the MetAmp directory and execute makefile: `make`. First, the installation script will check for required tools: R, Python and GCC. Then it builds and installs required packages and data for evaluation. You can also run `make install` to avoid checking stage.

3 Quick start

1. Edit script main.R, specifically:

- Provide your sequence library (amplicon sequence data), for example:

```
libs <- c("data/staggered/SRR072223_V13V31_1_relabeled.fastq",  
"data/staggered/SRR072223_V35V53_1_relabeled.fastq",  
"data/staggered/SRR072223_V69V96_1_relabeled.fastq")
```

- Provide reference sequences, for example:

```
# Reference 16S gene sequences
ref16S <- "data/16S_gold_hmc.fasta"
refs <- c("data/16S_gold_hmc_V13V31.fasta", # V1-3
"data/16S_gold_hmc_V35V53.fasta", # V3-5
"data/16S_gold_hmc_V69V96.fasta") # V6-9
```

2. Set the program directory in `main.R`:

```
dir_path <- " /Projects/metamp/"
```

3. Then run the script `main.R` from R-environment:

```
source("main.R")
```

3.1 Test data

The following test data sets were used:

- data/even/
- data/staggered/

The even and staggered are the Human mock community pyrosequence data (SRX021555), even and staggered community, 21 microbial species and one eucariotic.

3.2 Results

Results contain:

- Clustering results (`*.clstr`) - by default it is `clusters.clstr` but the name can be changed (by editing `clust_filename` parameter in `config.R` file)
- `otu_table.txt` - an OTU table, where rows are OTUs and columns are barcodes.
- A text file that contain coordinates for each point, including reference points (`coordinates.crd`)
- A text file containing all messages during the analysis process (`log.txt`)

4 MetAmp workflow

In the following sections we present key stages of our meta-amplicon analysis pipeline.

4.1 Input data

Input data can be any of the following NGS libraries: Roche 454, Ion Torrent, Illumina.

Roche 454 and Ion Torrent provide single-end sequences in special binary Sequence Flowgram Format (SFF). Before using SFF files you have to convert them into FASTQ format.

Illumina paired-end reads should be merged into longer single-end sequences.

4.2 Data denoising

Data denoising is necessary for downstream analysis because it removes the majority of 'foreign' nucleotides, such as barcodes and primers, and low-quality regions. We suggest you to use SeqyClean tool for this operation, since it is the most comprehensive tool for sequence denoising. You can download SeqyClean here: <https://bitbucket.org/izhbannikov/seqyclean>. Using SeqyClean is simple and straightforward. Nevertheless, we provided a special script that does all the data de-noising.

4.3 Analysis workflow

1. Building a reference topology of microbial populations, where pairwise distances are computed from applying pairwise alignment of complete reference 16S gene sequences from RDP database, see Figure 1(a), top plane. We use the percent sequence identity in order to compute the distances between sequences.
2. Computing a guide "amplicon" topology for each reference 16S sequence. To do this, we extract marker regions (for example V1- V3 and V6-V9) from reference complete 16S sequences ("reference points") and place them on to the plane with the same methodology that was used for building the reference topology of complete 16S (see Figure 1(a), bottom plane, hollow green circles).
3. add empirical amplicon sequences (in fact, amplicon consensus sequences), obtained through amplicon sequencing of microbial data (Figure 1b, bottom, filled blue circles). This topology is the same empirical topology as if the empirical data would be clustered with existing methods
4. Normalized empirical topologies are formed for each marker through mapping of each empirical point back to the reference 16S sequence topology. Mapping between the reference topology and guide "amplicon" topology is achieved with affine transformation. Guide sequences are so match with the corresponding reference 16S sequences, from where they were obtained from (see Figure 1a, top, solid green circles), and, in turn, carry the empirical amplicon sequences (or consensus amplicon sequences) back to the 16S plane.
5. Repeating stages 2-4 for each variable region.
6. Clustering and building OTUs.

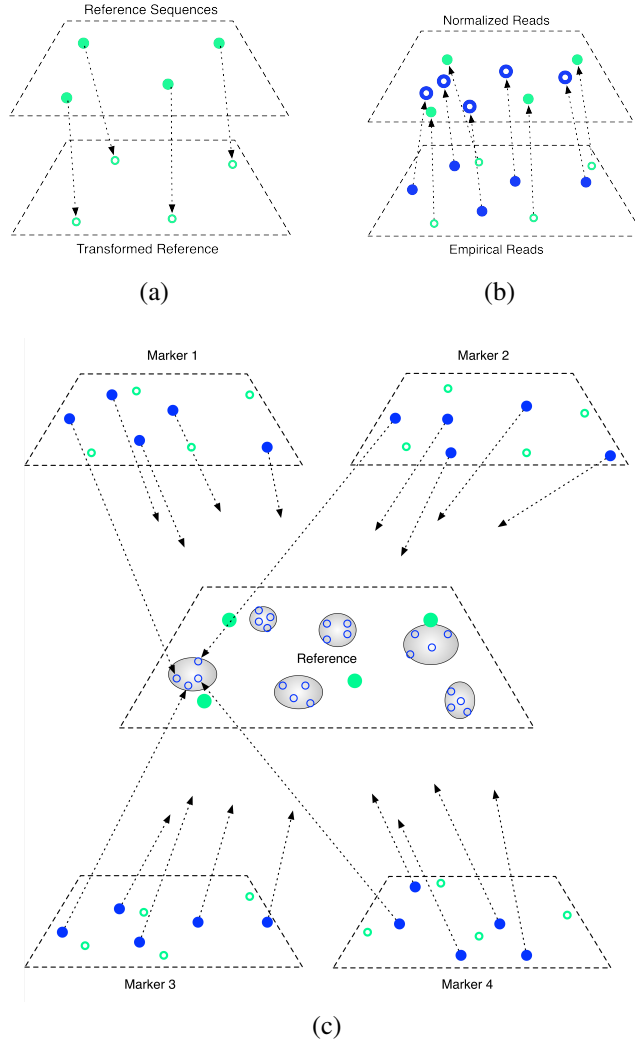


Figure 1: Illustration of “meta-amplicon” analysis algorithm. (a) Computing reference (top plane) and guide (bottom plane) topologies. Pairwise sequence dissimilarity is used to compute distances between reference sequences. The same approach is used for guide sequences, originated from corresponding whole 16S gene sequences. Reference and guide points are then placed on to the planes with multidimensional scaling. (b) The empirical amplicon sequences (hollow blue circles) are placed on to the bottom plane along with the guide circles (hollow green) circles and with the same metric. Such guide sequences are then mapped back to the reference 16S plane, carrying empirical sequences (solid blue circles). (c) Then this is repeated for each of marker (variable) region.

7. Final statistics.

4.4 Output files

The results contain:

- Clustering results (`*.clstr`) - by default it is `clusters.clstr` but the name can be changed (by editing `clust_filename` parameter in `config.R` file)
- `otu_table.txt` - an OTU table, where rows are OTUs and columns are barcodes.
- A text file that contain coordinates for each point, including reference points (`coordinates.crd`)
- A text file containing all messages during the analysis process (`log.txt`)

5 Provided data

We provided some data that can be used in your projects. Notice that the amount of data is very large (even unzipped files can take several Gb). For evaluation purposes we provided small data set. Description of data uploaded to our storage is given below:

- gold folder contains 'Gold Dataset' (`gold.fa`) containing the ChimeraSlayer reference database in the Broad Microbiome Utilities (<http://microbiomeutil.sourceforge.net/>), version `microbiomeutil-r20110519`. This contains all sequences in both orientations (forward and reverse-complemented.).
- LTP folder contains data (16S sequences, extracted marker regions and distance matrices) for data extracted from All Species Living Tree Project.

LTP:

| Name | Description |
|---------------------|------------------------------------|
| LTP-10271.fasta | Whole reference 16S sequences |
| LTP-10271_V13.fasta | reference sequences for marker 1-3 |
| LTP-10271_V35.fasta | sequences for marker 3-5 |
| LTP-10271_V69.fasta | sequences for marker 6-9 |

- Evaluation data set contains a small set of reference sequences (21 sequences) and empirical amplicon sequences (data/even.zip and data/staggered.zip). Amplicon sequences contain sequence data from Human Mock Community. More information provided here: <http://www.hmpdacc.org/HMMC/>

Evaluation:

| Name | Description |
|---------------------------|--|
| even | Even community: SRR053818, SRR072220, SRR072239 |
| staggered | Staggered community: SRR072221, SRR072223, SRR072237 |
| 16S_gold_hmc.fasta | contains 21 reference 16S sequences |
| 16S_gold_hmc_V13V31.fasta | contains 21 reference marker sequences (1-3 regions) |
| 16S_gold_hmc_V35V53.fasta | contains 21 reference marker sequences (3-5 regions) |
| 16S_gold_hmc_V69V96.fasta | contains 21 reference marker sequences (6-9 regions) |

The even and staggered are the Human mock community (HMP) pyrosequence data (SRX021555: <http://www.ncbi.nlm.nih.gov/sra?term=SRX021555>). Detailed description of these datasets (and sequencing protocols) is under the following link: <http://www.hmpdacc.org/HMMC/>


To run the program on test data open `main.R` and set the work directory (this directory should contain test data directory) and run the script `main.R` from R-environment:

```
>source("main.R")
```

Inside of the program directory we also provided some data.

5.0.1 Directory map and descriptions of folders

MetAmp

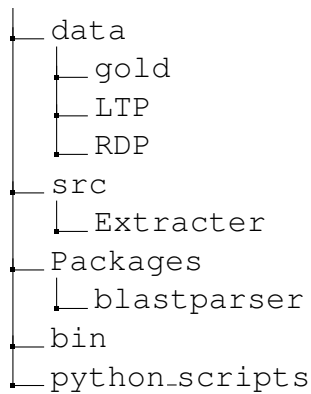

data - contains reference sequences, evaluation data and other data that can potentially be useful

gold - contains 'gold' reference sequences from Broad Institute. **data/LTP** - contains reference sequences and pre-computed distance matrices for > 10,000 species from All Species Living Tree Project.

src - contains source files.

Packages - contains installation files.

Packages/blastparser - a parser for blastn files.

bin - binary files (alignment and clustering programs).

python_scripts - Python scripts from www.drive5.com.

5.1 Hints on running your meta-amplicon analysis

Before analysis, you may need to perform re-labeling of your read headers: a read header should contain `barcodeLabel` (see more at www.drive5.com), for example:

```

@read1;barcodeLabel=TCAG;
GATGAACGCTGGCGGCGTGCCTAATACATGCAAGT...AT
+
IIIAAAIIIIIIIIIIIIIIIIHHHIAAAAAAAAAA...II

```

If you use Illumina sequence data, you also have to merge overlapping paired-end reads.

Later I will provide the scripts that can do these things above, but for now you can simply use scripts provided at www.drive5.com.

1. Edit script `main.R`, specifically:

Provide your sequences, for example:

```
libs <- c("data/staggered/SRR072223_V13V31_1_relabeled.fastq", # V1-3
          "data/staggered/SRR072223_V35V53_1_relabeled.fastq", # V3-5
          "data/staggered/SRR072223_V69V96_1_relabeled.fastq") # V6-9
```

and reference sequences:

```
# Reference 16S gene sequences
ref16S <- "data/16S_gold_hmc.fasta"
refs <- c("data/16S_gold_hmc_V13V31.fasta", # V1-3
          "data/16S_gold_hmc_V35V53.fasta", # V3-5
          "data/16S_gold_hmc_V69V96.fasta") # V6-9
```

2. Set the program directory in `main.R`, for example:

```
dir_path <- "~/Projects/metamp/"
```

3. Run the script `main.R`:

```
source("main.R")
```

6 Main parameters

Here we describe **config.R**:

analysis_dir - contains analysis files.

ref16S - path to the file with reference 16S sequences (file must be in FASTA format)

refs - path to the file with guide sequences, i.e. variable regions (file must be in FASTA format)

libs - data (empirical libraries)

7 Acknowledgements

This work was made possible by NIH Grants P20GM16448 and P20GM16448, and NSF Grant DBI0939454.