

September 14, 2014

MetAmp User Manual

Ilya Y. Zhbannikov¹, Janet E. Williams¹, James A. Foster^{1,2,3}

³Institute for Bioinformatics and Evolutionary Studies, University of Idaho, Moscow, ID 83844-3051,

²Department of Biological Sciences, University of Idaho, Moscow, ID 83844-3051,

¹Department of Bioinformatics and Computational Biology, University of Idaho, Moscow, ID 83844-3051,

Keywords:

Bioinformatics, metanomics, 16S genes, microbial studies, data mining

Correspondence:

Ilya Y. Zhbannikov

Department of Bioinformatics and Computational Biology

University of Idaho

Moscow, ID 83844-30521

Email: zhba3458@vandals.uidaho.edu

1 Description

MetAmp tool is developed for analysis of amplicon data by combining several marker regions from 16S rRNA genes. Such marker regions serve as unique identifiers for species. There are nine marker regions in total in bacterial 16S rRNA gene.

This tool developed by Ilya Y. Zhbannikov (zhba3458@vandals.uidaho.edu, i.zhbannikov@mail.ru), and James A. Foster (foster@uidaho.edu) with help of Janet Williams (janetw@uidaho.edu).

2 Installation

1. Download or clone MetAmp from the following repository:
<https://bitbucket.org/izhbannikov/metamp>
2. Execute 'Makefile' from MetAmp folder as follows: `$make`. First, the installation script will check for required software: R and GCC. Then it build and install required packages and data for evaluation. You can also run `$make install` to avoid checking stage.

3 Quick start

1. Edit script `config.R`, specifically: Provide your sequence library (amplicon data), for example:

```
libs <- c("<MetAmp directory>/test_data/Evaluation/Dataset1/SRR053817_V1_V3/SRR053817_3_cleaned.fasta",  
         "<MetAmp directory>/test_data/Evaluation/Dataset1/SRR053818_V3_V5/SRR053818_3_cleaned.fasta",  
         "<MetAmp directory>/test_data/Evaluation/Dataset1/SRR053819_V6_V9/SRR053819_3_cleaned.fasta")
```

and reference and guide sequences, for example:

```
ref16S <- "<MetAmp directory>/test_data/RDP/refSet16S-100.fasta" # Reference 16S gene sequences  
refs <- c("<MetAmp directory>/test_data/RDP/refSetV13-100.fasta", # Reference guide (marker) regions  
         "<MetAmp directory>/test_data/RDP/refSetV35-100.fasta",  
         "<MetAmp directory>/test_data/RDP/refSetV69-100.fasta")
```

If you do not provide them, the program will use previously aligned sequences:

```
res16S <- "<MetAmp directory>/test_data/RDP/res16S-100.usearch" # Pre-aligned reference 16S gene sequences.
refM <- c("<MetAmp directory>/test_data/RDP/resV13-100.usearch", # Pre-aligned guide (marker) regions.
        "<MetAmp directory>/test_data/RDP/resV35-100.usearch",
        "<MetAmp directory>/test_data/RDP/resV69-100.usearch")
```

2. Set the program directory in "main.R":

```
dir_path <- "~/Projects/metamp/"
```

3. Run the script "main.R":

```
source("main.R")
```

3.1 How to run the program on test data

The following test data sets were used:

- test_data/Evaluation/Dataset1.zip
- test_data/Evaluation/Dataset2.zip

The Dataset1.zip and Dataset2.zip are the Human mock community pyrosequence data (SRX021555), cleaned with SeqyClean.

To be able to use these sets you need to unzip them into any folder you want.

4 MetAmp workflow

In the following sections we present key stages of our meta-amplicon analysis pipeline.

4.1 Input data

Input data can be any of the following NGS libraries: Roche 454, Ion Torrent, Illumina.

Roche 454 and Ion Torrent provide single-end sequences in special binary Sequence Flowgram Format (SFF). Those reads are ready for using them in downstream analysis.

Illumina paired-end reads should be merged into longer single-end sequences. For this task we use FLASH program that provided by Magos and Salzberg (2011). For merging we provided an R script that performs this task.

4.2 Data denoising

Data denoising (also called "demultiplexing") is necessary for downstream analysis because it removes the majority of 'foreign' nucleotides, such as barcodes and primers, and low-quality regions. We suggest you to use SeqyClean tool for this operation, since it is the most comprehensive tool for sequence denoising. You can download SeqyClean here: <https://bitbucket.org/izhbannikov/seqyclean>. Using SeqyClean is simple and straightforward. Nevertheless, we provided a special script that does all the job.

4.3 Analysis workflow

1. Building a reference topology of microbial populations, where pairwise distances are computed from applying pairwise alignment of complete reference 16S gene sequences from RDP database, see Figure 1(a), top plane. We use the percent sequence identity in order to compute the distances between sequences.
2. Computing a guide "amplicon" topology for each reference 16S sequence. To do this, we extract marker regions (for example V1- V3 and V6-V9) from reference complete 16S sequences ("reference points") and place them on to the plane with the same methodology that was used for building the reference topology of complete 16S (see Figure 1(a), bottom plane, hollow green circles).
3. add empirical amplicon sequences (in fact, amplicon consensus sequences), obtained through amplicon sequencing of microbial data (Figure 1b, bottom, filled blue circles). This topology

is the same empirical topology as if the empirical data would be clustered with existing methods

4. Normalized empirical topologies are formed for each marker through mapping of each empirical point back to the reference 16S sequence topology. Mapping between the reference topology and guide "amplicon" topology is achieved with affine transformation. Guide sequences are so match with the corresponding reference 16S sequences, from where they were obtained from (see Figure 1a, top, solid green circles), and, in turn, carry the empirical amplicon sequences (or consensus amplicon sequences) back to the 16S plane.
5. Repeating stages 2-4 for each variable region.
6. Clustering and building OTUs.
7. Final statistics.

4.4 Output files

The results contain:

- Clustering results (*.clust) - by default it is clusters.clstr but the name can be changed (by editing final_clust_filename parameter in config.R file)
- A text file that contain coordinates for each point, including reference (coordinates.crd)
- A text file containing all messages during the analysis process (log.txt)

5 Provided data

We provided some data that can be used in your projects. We used data hosting provided by www.copy.com. Notice that the amount of data is very large (even unzipped files can take several

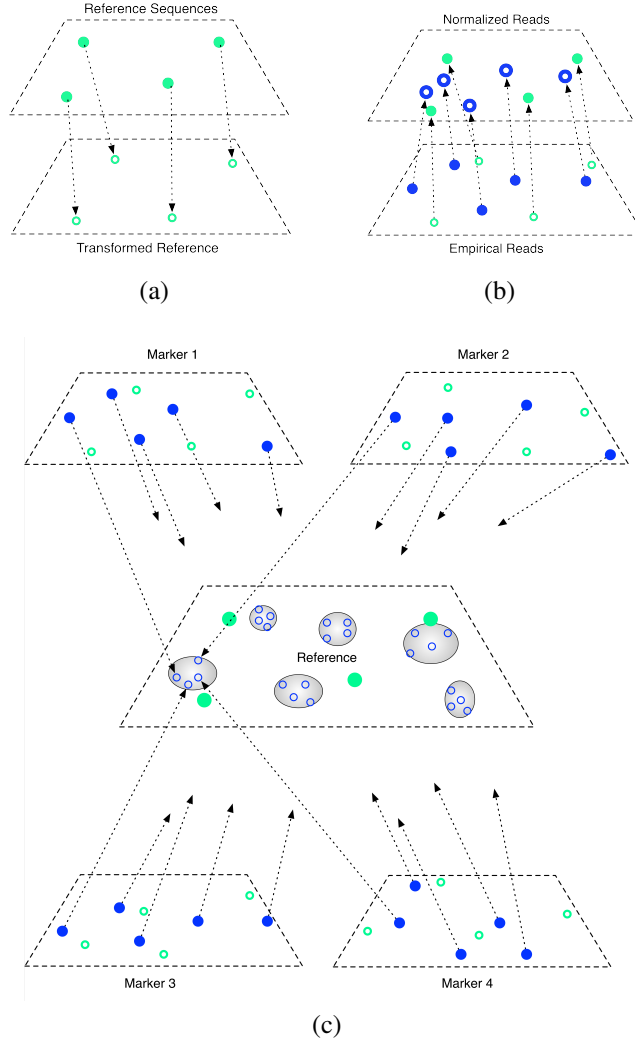


Figure 1: Illustration of “meta-amplicon” analysis algorithm. (a) Computing reference (top plane) and guide (bottom plane) topologies. Pairwise sequence dissimilarity is used to compute distances between reference sequences. The same approach is used for guide sequences, originated from corresponding whole 16S gene sequences. Reference and guide points are then placed on to the planes with multidimensional scaling. (b) The empirical amplicon sequences (hollow blue circles) are placed on to the bottom plane along with the guide circles (hollow green) circles and with the same metric. Such guide sequences are then mapped back to the reference 16S plane, carrying empirical sequences (solid blue circles). (c) Then this is repeated for each of marker (variable) region.

Gb). For evaluation purposes we provided small data set. Description of data uploaded to our storage is given below:

- LTP folder contains data (16S sequences, extracted marker regions and distance matrices) for data extracted from All Species Living Tree Project. The size of this folder is about 3.5 Gb. So be careful about your disc space.

LTP:

Name	Description
LTP-10271_V13.usearch.zip	zipped distance matrix for marker 1-3
LTP-10271_V35.usearch.zip	reference sequences for marker 3-5
LTP-10271_V69.usearch.zip	reference sequences for marker 6-9
LTP-10271.usearch.zip	reference sequences for the whole 16S gene
LTP-10271_V13.fasta	reference sequences for marker 1-3
LTP-10271_V35.fasta	sequences for marker 3-5
LTP-10271_V69.fasta	sequences for marker 6-9

Link to this folder: https://copy.com/ULDKi_jpPra6n3hvA

- RDP folder contains the same type of data but extracted from Ribosomal Database Project:
<http://rdp.cme.msu.edu/>

RDP:

Name	Description
RDP_matrix16S-1000.usearch	distance matrix for marker 6-9
RDP_matrix69-1000.usearch	distance matrix for marker 6-9
RDP_matrix35-1000.usearch	distance matrix for marker 3-5
RDP_matrix13-1000.usearch	distance matrix for marker 1-3
RFP_ref13-1000.fasta	reference sequences for marker 1-3
RDP_ref35-1000.fasta	reference sequences for marker 3-5
RDP_ref69-1000.fasta	reference sequences for marker 6-9
RDP_ref16S-1000.fasta	whole 16S reference sequences

RDP folder is much smaller than LTP and takes about 300 Mb unzipped. Link to this folder:

<https://copy.com/1rHaBXQsXM7SUtk2>

- Evaluation data set contains a small set of reference sequences (100 sequences), amplicon sequences and pre-clustered data sets. Amplicon sequences contain sequence data from Human Mock Community. More information provided here: <http://www.hmpdacc.org/HMMC/>

Evaluation:

Name	Description
Dataset1.zip	zipped sequence data for even community (ids)
Dataset2.zip	zipped sequence data for staggered community (ids)
HMC_ref_16S.fasta	contains 30 reference 16S sequences
HMC_ref_V13.fasta	contains 30 reference marker sequences (1-3 regions)
HMC_ref_V35.fasta	contains 30 reference marker sequences (3-5 regions)
HMC_ref_V69.fasta	contains 30 reference marker sequences (6-9 regions)
HMC_ref_16S.usearch	contains distance matrix for reference 16S sequences.
HMC_ref_V13.usearch	contains distance matrix for reference marker 1-3 sequences
HMC_ref_V35.usearch	contains distance matrix for reference marker 3-5 sequences
HMC_ref_V69.usearch	contains distance matrix for reference marker 6-9 sequences

Link to this folder: <https://copy.com/KZvTQvh01vgh8dAp>

The Dataset1.zip and Dataset2.zip are the Human mock community (HMP) pyrosequence data (SRX021555: <http://www.ncbi.nlm.nih.gov/sra?term=SRR053818>), cleaned with SeqyClean. Detailed description of these datasets (and sequencing protocols) is under the following link: <http://www.hmpdacc.org/HMMC/>

To run the program on test data open test.R and set the working directory (this directory should contain test_data directory) and run the script test.R from R-environment:

```
>source("test.R")
```

By default we provided pre-computed datasets that contain:

- Distance matrix (and actual reference sequences in FASTA format) computed out of reference 16S genes. These reference sequences were downloaded from All Species Living Tree Project web site(<http://www.arb-silva.de>). There are 10,271 sequences in total.

- Distance matrices from three (V1-3, V3-5 and V6-9) marker regions extracted from the reference set described above. So in total, there are 10,271 sequences for each marker region.

Inside of the program directory we also provided some data.

5.0.1 *Directory map and descriptions of folders*

```

MetAmp
├── Evaluation
│   ├── analysis
│   └── data
├── data
│   └── RDP
├── src
│   ├── Extracter .2 install
│   ├── blastparser
│   └── cdhitconverter
└── bin

```

Evaluation folder contains data for evaluation purposes.

Evaluation/analysis - pre-computed (clustered) sequences to speed-up evaluation.

Evaluation/data - contains sequence libraries for evaluation.

data - contains reference sequences and other data that can potentially be useful

data/RDP - contains reference sequences and pre-computed distance matrices for 100 species from Ribosomal Database Project.

src - contains source files.

install - contains installation files.

install/blastparser - a parser for blastn files.

install/cdhitconverter - converts CD-HIT's .clstr data.

bin - binary files (alignment and clustering programs).

5.1 Hints on running your meta-amplicon analysis

Before analysis, you may need to perform data de-noising and (if you use Illumina sequence data), merge overlapping reads. These have been already described in the section above.

Later I will provide the scripts that do it.

1. Edit script config.R, specifically:

Provide your sequences, for example:

```
libs <- c("~/Projects/metamp/test_data/Dataset1/SRR053817_V1_V3/SRR053817_3_cleaned.fasta", # V1-3
          "~/Projects/metamp/test_data/Dataset1/SRR053818_V3_V5/SRR053818_3_cleaned.fasta", # V3-5
          "~/Projects/metamp/test_data/Dataset1/SRR053819_V6_V9/SRR053819_3_cleaned.fasta") # V6-9
```

and reference sequences:

```
ref16S <- "~/Projects/metamp/test_data/1000/refSet16S-1000.fasta" # 16S gene sequences
refs <- c("~/Projects/metamp/test_data/1000/refSetV13-1000.fasta", # Marker regions
          "~/Projects/metamp/test_data/1000/refSetV35-1000.fasta",
          "~/Projects/metamp/test_data/1000/refSetV69-1000.fasta")
```

If you do not provide them, the program will use pre-aligned sequences:

```
res16S <- "1000/res16S-100.usearch" # Pre-aligned reference 16S gene sequences.
refM <- c("1000/resV13-100.usearch", # Pre-aligned guide (marker) regions.
          "1000/resV35-100.usearch",
          "1000/resV69-100.usearch")
```

2. Set the program directory in "main.R", for example:

```
dir_path <- "~/Projects/metamp/"
```

3. Run the script "main.R":

```
source("main.R")
```

6 Parameters

Here we describe **config.R**:

- analysis_dir** - contains analysis files.
- ref16S** - path to the file with reference 16S sequences (file must be in FASTA format)
- refs** - path to the file with guide sequences, i.e. variable regions (file must be in FASTA format)
- res16S** - pre-aligned reference 16S sequences in special 'uc' USEARCH format (see USEARCH manual for details)
- refM** - pre-aligned guide sequences in special 'uc' USEARCH format
- libs** - data (empirical libraries)
- denoise_app** - path to sequence denoising application (by default - SeqyClean)
- merge_app** - path to application that merges overlapping paired-end reads (by default - FLASH)
- use_custom_matrix** - F or T. 'T' indicates that the user wants to use pre-aligned references
- cluster_app** - clustering application (by default - CD-HIT). Clustering is used to reduce amount of data by computing consensus sequences.
- cluster_identity_threshold** - identity threshold for clustering application (by default - 0.9)
- cluster_suff** - suffix for clustered intermediate analysis files

7 Acknowledgements

This work was made possible by NIH Grants P20GM16448 and P20GM16448, and NSF Grant DBI0939454.